



THE HYGIENE SPECIALIST

A DIVERSE RANGE OF HYGIENE PRODUCTS

OUR PRODUCT SECTIONS

- HAND DISINFECTANT
- HAND DISPENSERS
- SKIN PREPARATION
- WASTE BINS
- PAPER DISPENSERS
- SURFACE & INSTRUMENTAL DISINFECTANT
- SURGICAL SCRUB SINK & TRANSFER TROLLEY

WWW.INOQUESTLABS.COM



CONTENTS

Company Profile	4
Hand Disinfectants	12
Skin Preparation	17
Surface & Instrument Disinfectants	24
Hand Dispensers	39
Surgical Sink & Transfer Trolley	49
Paper Dispensers	55
Waste Bins	64
Contact Us	70





COMPANY PROFILE

IQ is one of the leading distributors of Hospital
Care & Hygiene Products in United Arab Emirates

About Us

Established in Dubai 2010, Inoquest Labs Industries LLC (IQ) is committed to research, development and production of innovative products in disinfection, sanitation and cleaning.

It is our key concern to combine and balance our products efficiency with environmental responsibility. Inoquest Labs as a young and dynamic company - is always on the level of latest innovation and technological progress.

In addition to offering our products we are committed to develop with our customers a problem-solving relation. Our customers aim to have the best disinfectant and cleaning solutions are the focus of our efforts. With customized concepts and products we are partner to reach that goal on a daily basis. Through these unique efforts our ambition will be to build a durable trust and reliable cooperation with our customers.

66

Together with our customers,

For almost 30 years our management have been offering our business partners our expertise so that they can realise new product ideas in their name. Together with our customers, we develop product concepts that comply with the latest market requirements and legal regulations.



Our Vision

Position ourselves to be a regional player of Hospital Care & Hygiene Products. Improve on contemporary measures including cost, quality, service, speed of delivery and mobilization. Keep investing in technology, systems and human resource to effectively meet the challenges every new dawn brings.

Our mission

To be a partner in the markets of Health Care Industry by providing quality products and superior service to our customers, while learning from their feedback to set even higher standards for our products. IQ that endeavors to set the highest standards in corporate ethics in serving the society.

Our Value



Mutual respect binds us together. We comply with ethical work principles and encourage all of those we work with to follow suit. It matters to us how we nurture our relationships with employees, users and business partners. Every individual is important for IQ's success.



Together with our employees, we are building a safe, healthy and creative company. Our work is planned with people and the environment in mind. All resources are used prudently for a noble cause. By managing and constantly improving processes, we offer high-quality and innovative products.



We are everywhere, all of the time. Every touch is a unique experience. We invest in the development of people and equipment in cooperation with the best experts. We perform our work conscientiously and strive to be the best in everything we do.



Tablet Pressing Machine

Director Statements

IQ is a company that meets the hygiene needs of its customers day after day. Our hygienic products enhance quality of life. We offer comprehensive solutions to everyone who is a part of the IQ's story, and work with them to create opportunities for our mutual and sustainable development. Proudly since 2010.

Systematically deal with hospital micro-organisms
Multi-drug resistant micro-organisms are a known risk in hospitals and medical facilities. The risk of transmitting infection via personnel is high but can be prevented by effective hand disinfection. Regular surface disinfection and cleaning can also control the transmission of micro-organisms to patients and personnel.

We provide you with reliably effective and user-friendly products. For all medical and nursing areas.



Advanced Chemical Lab

Oliver Balk General Manager

Our many years of experience developing and manufacturing disinfectant and cleaning agents means that we provide you with reliably effective and user-friendly products for all medical and nursing areas.

As a professional partner, we support you with our extensive technical expertise, innovative products, a deep appreciation of quality and the necessary flexibility.

We are also happy to follow your specific requirements and turn your product ideas into reality.

Our Company

We create creative work to promote your business.

The success of any organization is mainly dependent on the people working for it. Each member of our team is considered equally important and provided constant training, motivation, and guidance. We possess a dedicated staff of the highest caliber committed to making our business a success.

Our policies and core values aim to create an informal yet stimulating team-oriented work environment with a culture of sharing and open communication. We cherish the diversity of viewpoint of every individual; we realize this encourages innovation and develops character. All employees have the right to a stress and injury-free work environment.

Our Products & Services

Our Product Portfolio is wide and will be continually extended. Markets are changing, and customer needs are changing. By implementing new products as per customer needs. We keep you in business. Products range from Hand Sanitizer, Medical and Hygiene Soaps, Shower Gel and Cosmetic, Instrument and Floor Cleaner, Dish Wash & Laundry Liquids, Detergents and Cleaning Agents, Dispensing Systems like wall Dispenser, etc., Scrub Sink & Waste Bin, Heel Incision Lancets & Medical Disposable.



**HOSPITAL & HYGIENE
PRODUCTS**



**DISINFECTANTS &
DETERGENTS**



**HOME CARE
PRODUCTS**



**DISINFECTION
EQUIPMENT**



**PRIVATE LABEL OR
OEM PRODUCTS**

Our Plant

IQ is a manufacturing company in Dubai with branches in Bahrain and Germany. We are specialized in Hygiene Solutions for Healthcare and other Hygiene-related industries.

IQ is also one of the leading distributors of Medical and General Waste Bins in United Arab Emirates and we are representing brands in the field of healthcare, hygiene and medical consumable like OPHARDT HYGIENE Germany, ORO- Chemie Germany, KANMED Sweden, MEDBAR Turkey and others.

Our Manufacturing plant:

DIN ISO 9001:2015 Certified Company

GMP & DIN ISO 22716:2007

7 Filling lines for liquids; Filling volume 50ml - 10 Ltr;


4 Mixer for liquids, Liquid production capacity per day 26,000 Liter.

Tablet Pressing Machine EURO B up to 25mm ; 20 Stations ;

Mixer for powder Production capacity of tablets: max. 120,000 per hour.

Storage Capacity: 600 pallets in high rack storage facility.



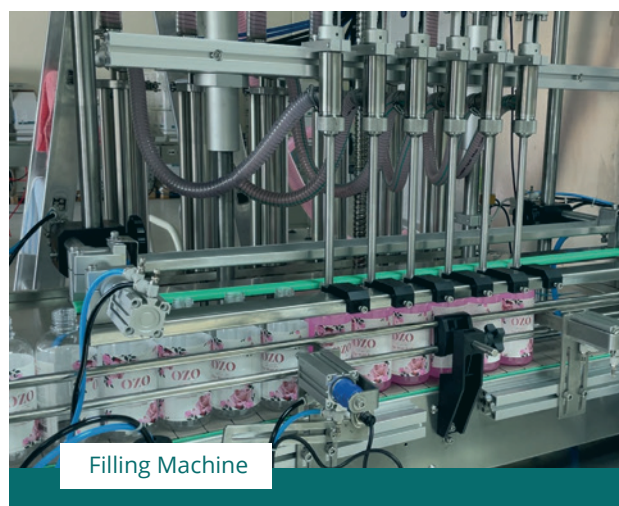



Loading Area



Mixing



Ready For Delivery



INOQUEST LABS
Disinfectant and Hygiene Products



HAND DISINFECTANTS

The majority of work in hospitals is performed by hand, meaning they also come into contact with countless bacteria, fungi, and viruses with the potential to cause diseases. The risk of transmitting infection via personnel is therefore very high but can be prevented with effective hand disinfection. Due to the frequency of skin disinfection in medical facilities, good compatibility of the products and regular care of the hands is important to maintain the health of the skin. IQ hand disinfectants have fast and reliable disinfection action while also offering very good skin compatibility.

IQ Hand Sanitizer

Scent Free – GEL

Isopropyl alcohol 70% v/v

The risk of transmitting infection via personnel in hospitals is very high but can be prevented with effective hand disinfection. IQ Hand Sanitizer Gel contains 70% of Ethanol alcohol and kills 99.99% of most common germs that may cause illness in as little as 30 seconds. Alcoholic gel to be rubbed into the skin. Precise dosing and easier, non-drip application through user-friendly gel formulation and dispensing systems. Use without water to supplement routine handwashing. Our Gel is free from cumulative long-term acting ingredients that might have an irritating or sensitizing effect on the skin. Compatible with washing lotion and latex gloves.

► IQ Pocket Hand Sanitizer

Easy to handle, easy to carry! Our pocket bottles are ideal to carry around for a fast and reliable source of proper hand disinfection.



Table Top	U/M ml	Pack Size	Art. No.
Pocket bottle with Flip Cap	50 ml	200 pcs/ box	I020480200
Alcoholic Gel Bottle with Flip Cap	150 ml	45 pcs/ box	I020860045
Alcoholic Gel Bottle with Flip Cap	200 ml	40 pcs/ box	I020860242

► IQ Hand Sanitizer Table Top

Our range for tabletop bottles is wide. From 250ml up to 1 Litre every bottle volume or design is possible. All bottles are coming with pumps.



Dispensers	U/M ml	Pack Size	Art. No.
Desktop Bottle with Pump	250 ml	14 pcs/ box	I020480014
Round Bottle with Pump	500 ml	06 pcs/ box	I020480016
Euro Bottle with Pump	1000 ml	08 pcs/ box	I020480083
Euro Bottle with Pump	500 ml	20 pcs/ box	I020480052
Square Bottle with Pump	500 ml	20 pcs/ box	I020480053

Active Ingredients: Aqua, Polymer, Ethanol 70%, Adjuvants.

► IQ Pouch Dispensers

Pouches are the most economical form of dispensing liquids in health care. We can provide IQ hand Disinfectant in different pouch systems. US and Chinese dispensers are the most common in the market.



Table Top	U/M ml	Pack Size	Art. No.
Round Nozzle Pouch Dispenser	1000 ml	12 pcs/ box	I020480112
Tubular Pouch Dispenser	1000 ml	12 pcs/ box	I020480012
Triangle Pouch Dispenser	1000 ml	10 pcs/ box	I020480110

► IQ Hand Sanitizer Hygienic Dispensers

We can provide IQ Hand Disinfectant nearly every common bottle for European, US, and Chinese wall dispensers. Economical with cap and as a cartridge system with pump.



Dispensers	U/M ml	Pack Size	Art. No.
Euro Bottle	1000 ml	08 pcs/ box	I020480088
Euro Bottle	500 ml	20 pcs/ box	I020480020
KX Bottle	1000 ml	12 pcs/ box	I020480212
Neptune Bottle	1000 ml	08 pcs/ box	I020480008

► IQ Hand Sanitizer 5L Can

Bulk quantities in can, drum, or IBC container for commercial use. Please ask for specific packaging solutions.



Table Top	U/M ml	Pack Size	Art. No.
IQ Hand Sanitizer 5L Can	5 Litre	04 pcs/ box	I020480054

Isopropyl alcohol 70% v/v



Liquid



Table Top



Dispensers

Active Ingredients : Aqua, Monopropylenglycol, Isopropyl alcohol 70% v/v.

IsoPropyl alcohol 70% v/v Chlorhexidine gluconate 0.5%



Liquid



Table Top



Dispensers

Active Ingredients:
Aqua, Glycerine, Monopropylenglycol, Chlorhexidine Gluconate, Isopropyl alcohol, Propanediol.

► IQ C20 Hand + Skin Disinfectant Hand + Skin Disinfectant

IQ C20 is a fast-acting, alcohol-based rub-in preparation for hygienic and surgical hand disinfection in clinics, care homes (old peoples homes, etc.), the food industry, and all fields impeccable hygiene conditions are necessary. IQ C20 is also suitable for disinfecting skin before simple injections and punctures of peripheral vessels, disinfecting the skin before operations and holes of joints, and disinfecting skin rich in sebaceous glands, as well as for cold compresses and for disinfecting suitable gloves.

Effective spectrum:

Bactericidal, tuberculocidal, fungicidal, virucidal against (enveloped viruses, e.g. Vaccinia viruses, HBV, HCV and HIV, SARS-Corona viruses as well as non-enveloped viruses, e.g. rotaviruses, noroviruses).

► IQ C20 Table Top

Our range for tabletop bottles is wide. From 120ml up to 1 Litre every bottle volume or design is possible. All bottles are coming with pumps.

Table Top	U/M ml	Pack Size	Art. No.
Euro Bottle with Pump	1000 ml	08 pcs/ box	1020090083
Euro Bottle with Pump	500 ml	20 pcs/ box	1020090052
Square Bottle with Pump	500 ml	20 pcs/ box	1020090053

► IQ C20 Hygienic Dispensers

We can provide IQ C20 in nearly every common bottle for European, US, and Chinese wall dispensers. Economical with cap and as a cartridge system with pump.

Dispensers	U/M ml	Pack Size	Art. No.
Euro Bottle	1000 ml	08 pcs/ box	1020090088
Euro Bottle	500 ml	20 pcs/ box	1020090020
KX Bottle	1000 ml	12 pcs/ box	1020090212
Neptune Bottle	1000 ml	08 pcs/ box	1020090008

► IQ ISO Hex 0.5% Antiseptic Hand Rub

IQ ISO Hex 0.5% is a skin antiseptic with a remanent effect. Used as a skin antiseptic before invasive/surgical operations and general hygiene applications. Effective against viruses, bacteria, and fungi. With its cationic structure, chlorhexidine gluconate bonds strongly with the tissue and continue to show its protective property for several hours after application. IQ ISO Hex 0.5% is well tolerated by the skin.

Kills 99.999% of bacteria, no perfume, no dye, latex glove compatible.

Antimicrobial Properties:

Virus inactivating (HBV, HCV, HIV), Bactericidal (incl. TB, MRSA, E-coli), Fungicidal IQ ISO Hex 0.5% contains isopropanol 70% and chlorhexidine gluconate 0.5% as active compound plus auxiliary substances.

► IQ ISO Hex 0.5% Table Top

Our range for table-top bottles is wide. From 120ml up to 1 Litre every bottle volume or design is possible. All bottles are coming with pumps.

Table Top	U/M ml	Pack Size	Art. No.
Euro Bottle with Pump	1000 ml	08 pcs/ box	1020650083
Euro Bottle with Pump	500 ml	20 pcs/ box	1020650052
Square Bottle with Pump	500 ml	20 pcs/ box	1020650053

► IQ ISO Hex 0.5% Hygienic Dispensers

We can provide IQ ISO Hex 0.5% in nearly every common bottle for European, US, and Chinese wall dispensers. Economical with cap and as a cartridge system with pump.

Dispensers	U/M ml	Pack Size	Art. No.
Euro Bottle	1000 ml	08 pcs/ box	1020650088
Euro Bottle	500 ml	20 pcs/ box	1020650020
KX Bottle	1000 ml	12 pcs/ box	1020650212
Neptune Bottle	1000 ml	08 pcs/ box	1020650008

Chlorhexidine gluconate 2% Ethanol 70%



Table Top



Dispensers

Active Ingredients : Aqua, Ethanol, Chlorhexidine Gluconate, Propanediol.

Chlorhexidine Gluconate 0.5% Ethanol 70%



Table Top



Dispensers

Active Ingredients : Aqua, Glycerine, Monopropylene Glycol, Chlorhexidine Gluconate, Ethanol, Propanediol.

► IQ HEXA Derm Advance

Skin Antiseptic

IQ HEXA Derm Advance is a skin antiseptic with a remanent effect. Used as a skin antiseptic before invasive/surgical operations and general hygiene applications. Effective against viruses, bacteria, and fungi. With its cationic structure, chlorhexidine gluconate bonds strongly with the tissue and continue to show its protective property up to 6 hr. after application. The absorption of chlorhexidine gluconate from the skin is almost negligible. IQ HEXA Derm Advance is well tolerated by the skin.

Antimicrobial Properties: Virus inactivating (HBV, HCV, HIV), Bactericidal (incl. TB, MRSA, E-coli), Fungicidal IQ HEXA Derm Advanced contains ethanol 70% and chlorhexidine gluconate 2% as active compound plus auxiliary substances.

► IQ HEXA Derm Advance Table Top

We can provide IQ HEXA Derm Advance in nearly every common bottle for European, US, and Chinese wall dispensers. Economical with cap and as a cartridge system with pump.

Table Top	U/M ml	Pack Size	Art. No.
Euro Bottle with Pump	1000 ml	08 pcs/ box	1020610083
Euro Bottle with Pump	500 ml	20 pcs/ box	1020610052
Square Bottle with Pump	500 ml	20 pcs/ box	1020610053

► IQ HEXA Derm Advance Hygienic Dispensers

Our range for tabletop bottles is wide. From 120ml up to 1 Litre every bottle volume or design is possible. All bottles are coming with pumps.

Dispensers	U/M ml	Pack Size	Art. No.
Euro Bottle	1000 ml	08 pcs/ box	1020610088
Euro Bottle	500 ml	20 pcs/ box	1020610020
KX Bottle	1000 ml	12 pcs/ box	1020610212
Neptune Bottle	1000 ml	08 pcs/ box	1020610008

► IQ HEXA Shield 0.5%

Antiseptic Hand Rub

IQ HEXA Shield 0.5% is an instant hand antiseptic with a remanent effect. Used as a skin antiseptic before invasive/surgical operations and general hygiene applications. Effective against viruses, bacteria, and fungi. With its cationic structure, chlorhexidine gluconate bonds strongly with the tissue and continue to show its protective property for several hours after application.

IQ HEXA Shield 0.5% is well tolerated by the skin. Kills 99.999% of bacteria, No perfume, No dye, Latex glove compatible

Antimicrobial Properties: Virus inactivating (HBV, HCV, HIV), Bactericidal (incl. TB, MRSA, E-coli), Fungicidal.

► IQ HEXA Shield 0.5% Table Top

Our range for tabletop bottles is wide. From 120ml up to 1 Litre every bottle volume or design is possible. All bottles are coming with pumps.

Table Top	U/M ml	Pack Size	Art. No.
Euro Bottle with Pump	1000 ml	08 pcs/ box	1020720083
Euro Bottle with Pump	500 ml	20 pcs/ box	1020720052
Square Bottle with Pump	500 ml	20 pcs/ box	1020720053

► IQ HEXA Shield 0.5% Hygienic Dispensers

We can provide IQ HEXA Shields 0.5% in nearly every common bottle for European, US, and Chinese wall dispensers. Economical with cap and as a cartridge system with pump.

Dispensers	U/M ml	Pack Size	Art. No.
Euro Bottle	1000 ml	08 pcs/ box	1020720088
Euro Bottle	500 ml	20 pcs/ box	1020720020
KX Bottle	1000 ml	12 pcs/ box	1020720212
Neptune Bottle	1000 ml	08 pcs/ box	1020720008

70% Alcohol Moisturizer



Table Top



Dispensers

Active Ingredients: Aqua, Polymer, Glycerine, Ethanol 70% v/v, Triethanolamine, Monopropylenglycol.

► Alcoholic Hand Disinfectant Gel

For Hygienic And Surgical Hand Disinfection

Alcoholic hand disinfectant gel for hygienic and Surgical hand disinfection to be rubbed into the skin.

- Bactericidal, tuberculocidal, yeasticidal
- Limited spectrum
- Virucidal activity
- With skin-care and
- Moisturizing agents
- Perfume-free
- No dye

► Alcoholic Hand Disinfectant Gel Table Top

Our range for tabletop bottles is wide. From 120ml up to 1 Litre every bottle volume or design is possible. All bottles are coming with pumps.

Table Top	U/M ml	Pack Size	Art. No.
Euro Bottle with Pump	1000 ml	08 pcs/ box	1020101083
Euro Bottle with Pump	500 ml	20 pcs/ box	1020101052
Square Bottle with Pump	500 ml	20 pcs/ box	1020101053

► Alcoholic Hand Disinfectant Gel Hygienic Dispensers

We can provide Alcoholic Gel in nearly every common bottle for European, US, and Chinese wall dispensers. Economical with cap and as a cartridge system with pump.

Dispensers	U/M ml	Pack Size	Art. No.
Euro Bottle	1000 ml	08 pcs/ box	1020101088
Euro Bottle	500 ml	20 pcs/ box	1020101020
KX Bottle	1000 ml	12 pcs/ box	1020101212
Neptune Bottle	1000 ml	08 pcs/ box	1020101008

85% Alcohol



Table Top



Dispensers

Active Ingredients: Aqua, Ethanol 85%v/v, Glycerine, Monopropylenglycol, Adjuvant.

► Alcoholic Hand Rub

For Hygienic And Surgical Hand Disinfection

Alcoholic hand rub for hand and skin disinfection to be rubbed into the skin

- Bactericidal
- Tuberculocidal
- Yeasticidal
- With skin-care and
- Moisturizing agents
- Perfume-free
- No dye

► Alcoholic Hand Rub Table Top

Our range for tabletop bottles is wide. From 120ml up to 1 Litre every bottle volume or design is possible. All bottles are coming with pumps.

Table Top	U/M ml	Pack Size	Art. No.
Euro Bottle with Pump	1000 ml	08 pcs/ box	1020770083
Euro Bottle with Pump	500 ml	20 pcs/ box	1020770052
Square Bottle with Pump	500 ml	20 pcs/ box	1020770053

► Alcoholic Hand Rub Hygienic Dispensers

We can provide Alcoholic Hand Rub in nearly every common bottle for European, US, and Chinese wall dispensers. Economical with cap and as a cartridge system with pump.

Dispensers	U/M ml	Pack Size	Art. No.
Euro Bottle	1000 ml	08 pcs/ box	1020770088
Euro Bottle	500 ml	20 pcs/ box	1020770020
KX Bottle	1000 ml	12 pcs/ box	1020770212
Neptune Bottle	1000 ml	08 pcs/ box	1020770008



INOQUEST LABS
Disinfectant and Hygiene Products



SKIN PREPARATION

Patient preoperative skin preparation can help you standardize practice across treatment areas and improve patient care.



4% Chlorhexidine
Gluconate

LONG LASTING PRE-OPERATIVE PROTECTION

IQ HEX WASH

4% Chlorhexidine Gluconate

Antimicrobial Hand and Skin Antiseptic

**Hand and skin antiseptic
with emollients.**

IQ HEX Wash cleans and disinfects your hand and skin with long lasting residual effect for pre-operative hand scrub.

- Residual effect for at least 6 hours¹.
- Chlorhexidine Gluconate with Isopropanol formula for added effectiveness.
- Slightly perfumed with emollients.
- IQ HEX WASH is also recommended for use as a pre-operative showering skin antiseptic prior to surgery.

What to know about HEX WASH?

LONG LASTING EFFECT

The combination of Chlorhexidine Gluconate with Isopropanol for added effectiveness for at least 6 hours. Chlorhexidine Gluconate builds a barrier on your skin keeps killing bacteria long after application.

BROAD SPECTRUM EFFICACY

The active chlorhexidine gluconate, a broad spectrum antiseptic, is effective against a wide range of both Gram positive and Gram negative bacteria (including MRSA), yeast, fungi and viruses².

SKIN FRIENDLY

With added emollients that keep the skin moist and protected, IQ HEX WASH can be used with confidence, ensuring healthcare professionals can fully comply with recommended scrub practices³.

GETTING YOUR SKIN READY FOR SURGERY IQ HEX WASH contains 4% of Chlorhexidine Gluconate (CHG). CHG is the recommended ingredient for pre-operative patient skin preparation via showering or bathing⁴.

LONG LASTING PRE-OPERATIVE PROTECTION.

INDICATIONS FOR USE

IQ HEX WASH is an antimicrobial hand and skin cleansing antiseptic. It is used as a pre-operative hand scrub and pre-operative skin cleanser prior to surgery. Suitable for soiled hand and skin. Use with water.

INGREDIENTS

Aqua, Cocamidopropyl betaine, C19-15 Pareth-8 surfactants, Emollient, 4% Chlorhexidine Gluconate, Isopropanol, Color, Fragrance.

REFERENCES

1. Larson E (1995 1 APIC Guidelines.
2. Denton, GW; 'Chlorhexidine.' Chapter 15 in Disinfection. Sterilization and Preservation. Ed. Block SS. Fifth Ed. Lippincott Williams and Wilkins. 2001.
3. Association for Perioperative Practice (AIPP) A Guide to Surgical Hand Antisepsis October 2014.
4. J.E. Coia, G.J. Duckworth, D.I. Edwards, M.Farrington, C. Fry, H. Humphreys, C. Mallaghan, D.R. Tucker for the Joint Working Party of the British Society of Antimicrobial Chemotherapy, the Hospital Infection Society, and the Infection Control Nurses Association Guidelines for the control and prevention of meticillin-resistant Staphylococcus aureus (MRSA) in healthcare facilities. Journal of Hospital Infection 2006, 63S:S1—S44. Available at www.his.org.uk/itil_documents/MRSA_Guidelines_PDF.pdf (accessed 28 March 2007).

► HEX WASH Hygienic Dispensers

We can provide IQ HEX WASH in nearly every common bottle for European, US and Chinese wall dispenser. Economical with cap and as a cartridge system with pump.

Name	U/M ml	Pack Size	Art. No.
Euro Bottle	1000 ml	08 pcs/ box	I020420088
Euro Bottle	500 ml	20 pcs/ box	I020420083
KX Bottle	1000 ml	12 pcs/ box	I020420212
Neptune Bottle	1000 ml	08 pcs/ box	I020420008



► HEX WASH Table Top

Our sortiment for table top bottles is wide. From 120ml up to 1 Litre every bottle volume or design is possible. All bottles are coming with pumps.

Name	U/M ml	Pack Size	Art. No.
Euro Bottle with Pump	1000 ml	08 pcs/ box	I020420083
Euro Bottle with Pump	500 ml	20 pcs/ box	I020420052
Square Bottle with Long Nozzle Pump	500 ml	20 pcs/ box	I020420070
Small Bottle	120 ml	32 pcs/ box	I020420032



► HEX WASH Pouch

Pouches are the most economical form of dispensing liquids in health care. We can provide IQ HEX WASH in different pouch system. US and Chinese dispensers are the most common in the market.

Name	U/M ml	Pack Size	Art. No.
Round Nozzle Dispenser	1000 ml	12 pcs/ box	I020420112



► HEX WASH Sponge and Scrub Brush

Our approved formula of IQ HEX WASH is also available as a saturated sponge for direct use with water or as a pre-saturated scrub brush with nail pick.

Name	U/M ml	Pack Size	Art. No.
HEX WASH Sponge 4%	1 Pc	30 pcs/ box	I020420088
Chlorhexidine Scrub Brush 4%	1 Pc	40 pcs/ box	I020420083





7.5% Povidone-Iodine
USP

IQ BETA SCRUB

ANTISEPTIC CLEANSER FOR PREOPERATIVE
WASHING AND SCRUBBING

Active ingredient: Povidone-Iodine USP 7.5%

IQ Beta Scrub contains povidone iodine 7.5g per 100ml with a content of 10% available iodine as active compound plus auxiliary substances. For hand and skin wherever a topical anti microbial scrub is required on the hospital ward.

Applications:

Ready-to-use povidone-iodine solution
Rapid activity within 30 seconds against germs
Can be used undiluted or diluted for irrigation
Active against bacteria incl. MRSA, mycobacteria, fungi, enveloped viruses (HBV / HIV), protozoa, etc

Advantages:

For painless disinfection
Immediate disinfecting action and good residual effect
Creamy foam for a fine lather with small bubbles
Broad spectrum activity

For external use only.

When using this product do not use near the eyes. In case of contact, rinse eyes thoroughly with plenty of water.

Store in a cool dry place.

Do not store above 25°C.

Keep away from children.

Ingredients:

Active ingredient: Povidone-Iodine USP 7.5%

► IQ Beta Scrub Table Top

Our sortiment for table top bottles is wide. From 120ml up to 1 Litre every bottle volume or design is possible. All bottles are coming with pumps.



Name	U/M ml	Pack Size	Art. No.
Euro Bottle with Pump	1000 ml	08 pcs/ box	I020520083
Euro Bottle with Pump	500 ml	20 pcs/ box	I020520052
Square Bottle with Long Nozzle Pump	500 ml	20 pcs/ box	I020520070
Small Bottle	120 ml	32 pcs/ box	I020520032

► IQ Beta Scrub Dispensers

We can provide IQ Beta Scrub in nearly every common bottle for European, US and Chinese wall dispenser. Economical with cap and as a cartridge system with pump.

Name	U/M ml	Pack Size	Art. No.
Euro Bottle	1000 ml	08 pcs/ box	I020520088
Euro Bottle	500 ml	20 pcs/ box	I020520083
KX Bottle	1000 ml	12 pcs/ box	I020520212
Neptune Bottle	1000 ml	08 pcs/ box	I020520008





2% Chlorhexidine
Gluconate

IQ HEXA DUO

ANTIMICROBIAL HAND WASH

IQ HEXA Duo contains chlorhexidine gluconate 2% as active compound.

IQ HEXA Duo is an antimicrobial hand wash solution. With 2% Chlorhexidine Gluconate, it has a residual effect and remains active after rinsing and drying of hands. IQ HEXA Duo contains antimicrobial, skin-care and moisturising agents for best protection of your hands. IQ HEXA Duo is bacteriostatic, virostatic (enveloped viruses) and fungistatic.

Application:

Wet skin with plenty of water. Apply the required amount of soap, wash gently and rinse thoroughly with water. If necessary use a brush and or repeat procedure.

For external use only.

When using this product do not use near the eyes. In case of contact, rinse eyes thoroughly with plenty of water.

Store in a cool dry place.

Do not store above 25°C.
Keep away from children.

Ingredients:

Aqua, 2% Chlorhexidine Gluconate, Non-ionic and Cationic Surfactants, Isopropanol, Emollient

Allergic or anaphylaxis reaction can occur in response to Chlomexidine Gluconate.

► IQ Hexa Duo Table Top

Our sortiment for table top bottles is wide. From 120ml up to 1 Litre every bottle volume or design is possible. All bottles are coming with pumps.



Name	U/M ml	Pack Size	Art. No.
Euro Bottle with Pump	1000 ml	08 pcs/ box	I020540083
Euro Bottle with Pump	500 ml	20 pcs/ box	I020540052
Square Bottle with Long Nozzle Pump	500 ml	20 pcs/ box	I020540070
Small Bottle	120 ml	32 pcs/ box	I020540032

► IQ Hexa Duo Dispensers

We can provide IQ Hexa Duo in nearly every common bottle for European, US and Chinese wall dispenser. Economical with cap and as a cartridge system with pump.

Name	U/M ml	Pack Size	Art. No.
Euro Bottle	1000 ml	08 pcs/ box	I020540088
Euro Bottle	500 ml	20 pcs/ box	I020540083
KX Bottle	1000 ml	12 pcs/ box	I020540212
Neptune Bottle	1000 ml	08 pcs/ box	I020540008





10% Povidone-Iodine USP
Without Detergent

IQ BETA STANDARD

ACTIVE AGAINST BACTERIA, VIRUSES, PROTOZOA AND FUNGI

Active ingredient: 10% Povidone-Iodine (1% available iodine)

Powerful Antiseptic and Disinfectant.

IQ Beta Standard is powerful antiseptic and disinfectant. In contact with skin and mucous membranes the povidone iodine complex gradually releases the active iodine, which works quickly and has a microbicidal effect (effective against a broad spectrum of bacteria, viruses, fungi and protozoa).

- Preoperative and post-operative skin antiseptics.
- Preoperative antiseptic preparation of the lining of the vagina.
- Antiseptic preparation of the skin and mucous membranes of injections, aspiration needle biopsy.
- Maintaining oral hygiene before, during and after dental or oral surgery.
- Disinfect skin (prevention of bacterial, viral and fungal infections and protozoal infections).
- Antiseptic treatment of wounds (lacerations, abrasions and static decubital infected ulcers, pyoderma, secondarily infected dermatoses, burns of second and third degree).

Do not use IQ Beta Standard

- If you are hypersensitive (allergic) to iodine;
- If you have a dysfunction of the thyroid gland (eg nodular colloid goitre, endemic goitre and thyroiditis lymphocitic)
- If you have renal failure (renal failure) and or hepatic (liver failure)
- The infant body weight below 1500 g.

Be careful when applying

IQ Beta Standard solution is intended for external use. If the patient have damage to the kidney and or liver function, be special caution in frequent use and only strictly prescribed by the doctor.

Long-term and frequent use should be avoided in patients with burns, especially those that take up more than 20 % of the body surface as well as in patients with large open wounds. In the case of local irritation, redness and swelling of the solution application to the skin should be stopped. Frequent use of the solution to the skin should be avoided if you are taking medication at the same time containing lithium.

Long-term and frequent use should be avoided in patients with burns, especially those that take up more than 20 % of the body surface as well as in patients with large open wounds. In the case of local irritation, redness and swelling of the solution application to the skin should be stopped. Frequent use of the solution to the skin should be avoided if you are taking medication at the same time containing lithium.

It is noted that Povidone iodine 10 % solution can have spermicidal (destroying sperm) effect, and can not be applied to mucous membranes of the vagina if the patient wants to conceive. If changes in the oral mucosa, despite treatment does not heal for a period of 14 days again doctor must be informed. The product is not recommended to apply after animal bites or insect bites.

Interactions

An iodine absorption through the skin can affect the results of functional tests for the thyroid gland. Inadvertent contamination for skin tests for occult blood (blood traces in the stool), or blood in the urine may result in false positive findings. Co-administration of a solution of Silver sulfadiazine or chlorhexidine may lead to partial inactivation of the medicine.

PREGNANCY AND LACTATION

During pregnancy and lactation it should only be used if strictly prescribed by the doctor. This applies also to newborns and infants. Iodine absorbed from the solution to the skin can pass the placental barrier and is excreted in breast milk. Frequent use of the solution to the skin during pregnancy should be avoided because of the potential risk to the fetus (development and function of the thyroid gland of the fetus). If necessary Povidone iodine 10 % solution applied to the skin for a long period of time, nursing mothers are advised to stop breast feeding and switch to artificial infant feeding.

ACTIVE INGREDIENTS

10 g Povidone-Iodine USP per 100 ml solution (1% Iodine)

► IQ Beta Standard Table Top

Our sortiment for table top bottles is wide. From 120ml up to 1 Litre every bottle volume or design is possible. All bottles are coming with pumps.



Name	U/M ml	Pack Size	Art. No.
Euro Bottle with Pump	1000 ml	08 pcs/ box	I020640083
Euro Bottle with Pump	500 ml	20 pcs/ box	I020640052
Square Bottle with Long Nozzle Pump	500 ml	20 pcs/ box	I020640070
Small Bottle	120 ml	32 pcs/ box	I020640032

► IQ Beta Standard Dispensers

We can provide IQ Beta Standard in nearly every common bottle for European, US and Chinese wall dispenser. Economical with cap and as a cartridge system with pump.



Name	U/M ml	Pack Size	Art. No.
Euro Bottle	1000 ml	08 pcs/ box	I020640088
Euro Bottle	500 ml	20 pcs/ box	I020640083
KX Bottle	1000 ml	12 pcs/ box	I020640212
Neptune Bottle	1000 ml	08 pcs/ box	I020640008



INOQUEST LABS
Disinfectant and Hygiene Products



SURFACE & INSTRUMENT DISINFECTANTS

☎ 8696 450 4 00971

✉ info@inoquestlabs.com

🌐 www.inoquestlabs.com

► IQ ALCO



IQ ALCO SPRAY

Ready-to-use alcoholic solution for disinfection of surfaces

Adehyde-free ready-to-use solution for rapid disinfection as well as for cleaning & disinfection of sensitive surfaces of medical devices (operating tables, medical devices, dressing trolleys, stethoscopes, etc.) and other surfaces, particularly in poorly accessible areas, in clinics, in care homes, the food industry and in all fields where impeccable hygienic conditions are necessary.



IQ ALCO Mist Sprayer
Ready-to-use 1L



IQ ALCO Spray
Ready-to-use 1L

Microbiological effectiveness

IQ Alco has the following characteristics of performance:

Application	Time
General Surface Disinfection	1 min.
Bacteria inclusive of MRSA ₁ and fungi ₂	1 min.
TB bacteria	1 min
Vaccinia viruses, BVDV incl. HBV, HCV and HIV ₃ , H1Ni, H5N2	2 min

IQ ALCO WIPES

Rapid alcohol disinfectant wipes for cleaning and disinfection of surfaces



IQ ALCO Wipes
100 Sheets



IQ ALCO Wipes
200 Sheets

Aldehyde-free concentrate for simultaneously disinfecting and cleaning all washable surfaces and objects (fixtures, devices) in all fields in clinics, e.g. in operating rooms, intensive care units, casualty wards, central sterilization departments and for diagnosis devices which are also used for operations.

IQ ALCO Wipes is also used by rescue services, in care homes (old people's homes, etc.), and in all fields where impeccable hygienic conditions are necessary

Microbiological effectiveness

IQ ALCO Wipes has the following characteristics of performance:

Application	Time
Surface disinfection with high load	5 min.
Bacteria inclusive	1 min.
Rota Virus	1 min
Papova & Polyoma viruses	2 min
Vaccinia viruses, BVDV incl. HBV, HCV and HIV ₃ , H1Ni, H5N2	1 min
Noro viruses	2 min

Product Name	U/M ml/sheets	Pack Size	Art. No.
ALCO with mist Sprayer	500 ml	20 pcs/ box	1020730056
ALCO Spray	1000 ml	09 pcs/ box	1020730039
ALCO Wipes	100 Sheets	20 pcs/ box	1070090012
ALCO Wipes	200 Sheets	12 pcs/ box	1070090020

► IQ B20

IQ B20 WIPES

Alcohol-free disinfecting wipes for cleaning and disinfection of alcohol-sensitive surfaces

Alcohol & aldehyde-free ready to use solution for simultaneously disinfecting and cleaning of alcohol-sensitive surfaces and devices and for all washable surfaces and objects in hospitals, medical and dental practices and ambulances.

IQ B20 ready-to-use solution is also used on surfaces in medical laboratories, by rescue services, in care homes (old people's homes, etc.), the food industry and in all fields where impeccable hygienic conditions are necessary.

Microbiological effectiveness

IQ B20 Wipes has the following characteristics of performance:

Application	Time
Surface disinfection with high load	5 min.
Bacteria inclusive	1 min
Rotavirus	1 min
Papova & Polyoma viruses	2 min
Vaccinia viruses, BVDV incl. HBV, HCV and HIV ₃ , H1Ni, H5N2	1 min



IQ B20 Wipes
100 Sheets



IQ B20 Wipes
200 Sheets

IQ B20 SPRAY

Alcohol-free disinfecting wipes for cleaning and disinfection of alcohol-sensitive surfaces

Alcohol & aldehyde-free ready to use solution for simultaneous disinfecting and cleaning of alcohol-sensitive surfaces and devices and for all washable surfaces and objects in hospitals, medical and dental practices and ambulances. IQ B20 ready-to-use solution is also used on surfaces in medical laboratories, by rescue services, in care homes (old people's homes, etc.), the food industry and in all fields where impeccable hygienic conditions are necessary.



**B20
Concentrate 1L**



**B20 Ready-to-use
Mist Sprayer 500ml**



**B20 Ready-to-use
Spray 1L**



**B20 Concentrate
Jerry Can 5L**

Microbiological effectiveness

IQ B20 Wipes has the following characteristics of performance:

Application	Time
Surface disinfection with high load	5 min.
Bacteria inclusive	1 min
Rotavirus	1 min
Papova & Polyoma viruses	2 min
Vaccinia viruses, BVDV incl. HBV, HCV and HIV ₃ , H1Ni, H5N2	1 min
Noro viruses	30 min



**B20 Concentrate
Stackable Jerry Can 10L**

Product Name	U/M ml/Sheets	Pack Size	Art. No.
B20 Wipes	100 Sheets	20 pcs/ box	I020090088
B20 Wipes	200 Sheets	12 pcs/ box	I020090020
B20 Concentrate	1000 ml	06 pcs/ box	I020100056
B20 Mist Sprayer	500 ml	20 pcs/ box	I020100039
B20 Spray	1000 ml	09 pcs/ box	I020100006
B20 Concentrate Jerry can 5L	5000 ml	04 pcs/ box	I020100054
B20 Concentrate Stackable Jerry Can 10L	10000 ml	Single	I020100010

► IQ IPA

IQ IPA SPRAY

70% Isopropanol Disinfectant Solution
for cleaning and disinfection of surfaces

IPA Spray ready-to-use solution for rapid and effective disinfection of all kinds of alcohol resistant inventory and surfaces.

Isopropanol Spray ready-to-use solution for rapid and effective disinfection of all kinds of alcohol resistant inventory and surfaces. Particularly suitable in poorly accessible areas and difficult to reach surfaces in all fields where impeccable hygienic conditions are necessary. IQ IPA Spray is free from perfume and residue-free.



**IPA Ready-to-use
Spray Bottle 1L**



**IPA Ready-to-use
Spray Bottle 1L**

Microbiological effectiveness

IQ IPA has the following characteristics of performance:

Application	Time
General Surface Disinfection (clean condition)	1 min.
Bacteria inclusive (clean surface) Staphylococcus aureus, Escherichia coli, Pseudomonas aeruginosa & Enterococcus hirae.	1 min

IQ IPA WIPES

70% v/v alcohol disinfectant wipes cleaning



IPA Wipes 100 sheets



IPA Wipes 200 Sheets

Aldehyde-free disinfectant wipes, ready to use, for all washable surfaces and objects (fixtures, devices) in all fields in clinics, e.g. in operating rooms, intensive care units, casualty wards, central sterilization departments and for diagnosis devices. All surfaces and objects must be alcohol resistant.

IQ IPA Wipes is also used by rescue services, in care homes (old people's homes, etc.), and in all fields where impeccable hygienic conditions are necessary

Microbiological effectiveness

IQ IPA Wipes has the following characteristics of performance:

Application	Time
General Surface Disinfection	1 min.
Bacteria and Fungi	30 sec
MRSA	30 sec
TB Bacteria	30 sec
Vaccinia viruses, BVDV incl. HBV, HCV and HIV ₃ , H1Ni, H5N2	15 sec
Noro viruses	5 min
Rota viruses	30 sec
SARS Corona viruses	30 sec

Product Name	U/M ml/Sheets	Pack Size	Art. No.
IPA Spray	1000 ml	09 pcs/ box	I020210039
IPA 5L Can	5000 ml	04 pcs/ box	I020210054
IPA Wipes	100 Sheets	20 pcs/ box	I070070012
IPA Wipes	200 Sheets	12 pcs/ box	I070070020



► IQ A20

IQ A20 SPRAY

Preliminary Disinfectant and Cleaning of Instrument

IQ A20 is an aldehyde-free concentrate for simultaneously disinfecting and cleaning the general and surgical instruments made of metal, plastic materials, rubber and glass in clinics, care homes (old people's homes, etc.) or medical research, etc.

Advantages of the product

- Manifold fields of application – suitable for all instruments, including rotary instruments which are sensitive to alkali and alcohol.
- Tested effectiveness with comprehensive effective spectrum against bacteria, fungi and all viruses (virucidal against enveloped and non-enveloped viruses).
- VAH/DGHM listed.
- Tested in accordance with EN 13727, EN 13624, EN 14348, EN 14476, EN 14561, EN 14562, EN 14563.
- Economical due to low application concentration.
- Extraordinary cleaning activity – removes the most stubborn deposits.
- Special corrosion inhibitors result in an excellent material compatibility.
- Toxicologically safe – the product does not contain any poisonous substances.
- Service life of ready-to-use solution: 14 days max.
- Pleasant, fresh fragrance.

Concentrate: 4 Years

Ready to use solution: * Fresh solution 28 days

* Used solutions 14 days max.

(Depending on the contamination)

Physical Properties

Appearance		Density	pH value
Concentrate	Clear green liquid of low viscosity	$D = 1.05 \pm 0.01 \text{ g/cm}^3 (20^\circ\text{C})$	13.8 ± 0.5
1% Ready-to-Use Solution	Clear, light green solution	—	11.9 ± 0.5

Toxicology

The acute oral toxicity (LD50) of the concentrate in rats is approx. 1400 mg/kg of body weight. IQ A20 is therefore classified as "harmful". The acute dermal toxicity (LD50) in rats is > 2000 mg/kg of body weight. IQ A20 is therefore classified as "non-toxic". Contact of the 2% ready-to-use solution with the eyes causes irritation whereas contact with the skin does not have any irritating effects.

Ecology

If adequately diluted, the ready-to-use solutions are ecologically safe as all organic constituents are biodegradable. The packing is made of polyethylene and can therefore be exploited both materially and thermally. To recycle, rinse the containers.

Microbiological effectiveness

IQ A20 Instrument Disinfectant has the following characteristics of performance:



**A20 Ready-to-use
Spray 1L**

Application	Concentration	Time
Instrument disinfection (according to VAH/DGHM)	2%	15 min.
Bacteria inclusive of MRSA ₁ and fungi ₂	2%	15 min.
TB bacteria	2%	1 hour
Vaccinia viruses, BVDV incl. HBV, HCV and HIV ₃	1%	1 hour
Adeno viruses	2%	30 min.
Polyoma SV 40 viruses	2%	5 min.
Polio viruses	6%	4 hour
Noro viruses ₄	1%	4 hour

- 1) Effectiveness against bacteria implies effectiveness against MRSA.
- 2) Test for *Candida albicans* and *Aspergillus niger*.
- 3) In accordance with RKI recommendation (Bundesgesundheitsbl. 47, 62 – 66, 2004).
- 4) Effectiveness against Polio viruses implies effectiveness against Noro viruses.



A20 Concentrate 1L



**A20 Concentrate Jerry can
Can 5L**



**A20 Concentrate Stackable
Jerry can 10 L**

Product Name	U/M ml	Pack Size	Art. No.
A20 Ready-to-use Spray	1000 ml	09 pcs/ box	I020060039
A20 Concentrate 1L	1000 ml	06 pcs/ box	I020060006
A20 Concentrate Jerry Can 5L	5000 ml	04 pcs/ box	I020060054
A20 Concentrate Stackable Jerry Can 10L	10000 ml	Single can	I020060010

► IQ ENZYMATIC CLEANER



IQ ENZYMATIC CLEANER

Concentrate for the chemothermic and manual reprocessing of instruments

IQ Enzymatic Cleaner is an enzymatic ready-to-use spray for the chemothermic and manual reprocessing of endoscopes and medical instruments. The product is characterized by very good soil-dissolving properties, which ensures a thorough cleaning. At the same time it is gentle on the material of endoscopes and medical instruments due to its neutral pH-value. IQ Enzymatic Cleaner is suitable for the use in hospitals, medical practices and laboratories and is also suitable for the use in ultrasonic baths (35°C - 45°C).

- *Very effective and efficient cleaning*
- *Free of perfumes and colourants*
- *Residue free in 2 - 3 minutes*
- *Good material compatibility*
- *Low foaming*
- *Bio- degradable*
- *Very economic use*
- *pH - neutral formula*

Application:

Manual reprocessing:

Dilute concentrate with water; choose a concentration between 0.5 and 3% depending on the degree of soiling. Completely immerse dismantled and opened instruments fully into the prepared solution for 2 - 3 minutes. Ensure that the entire instrument and parts are fully submerged during the entire contact time and avoid bubbles and cavities. Cover cleaning trays. The recommended water temperature is between 30°C and 40°C. Renew solution once or twice daily; if heavily contaminated, prepare fresh solution immediately.

Chemothermic reprocessing:

The dosage of 0.5 - 1.0% is regulated via the pumps inside endoscope cleaning and disinfection machines or other reprocessing machines. The recommended temperature for the water is between 35 C and 45 C. Follow the recommendations given by the manufacturer of the machines. After the cleaning step properly rinse instruments and parts. Instrument disinfecting procedure can be applied after rinsing cleaned instruments and parts.



**Enzymatic Cleaner
Ready-to-use
Spray Bottle 1L**



**Enzymatic Cleaner
Concentrate 1L**



**Enzymatic Cleaner
Concentrate Can 5L**

Product Name	U/M ml	Pack Size	Art. No.
Enzymatic Cleaner Ready-to-use	1000 ml	09 pcs/ box	1020110039
Enzymatic Cleaner Concentrate Round bottle	1000 ml	06 pcs/ box	1020110006
Enzymatic Cleaner Concentrate Jerry	5000 ml	04 pcs/ box	1020110054



► IQ B45

IQ B45 SPRAY

Quick - Acting Disinfectant
Adehyde-free ready-to-use solution

Adehyde-free ready-to-use solution for fast-acting alcoholic wipe or spray disinfection as well as for cleaning sensitive surfaces of medical devices (operating tables, medical devices, dressing trolleys, stethoscopes, etc.) and other surfaces, particularly in poorly accessible areas, in clinics, in care homes, the food industry and in all fields where impeccable hygienic conditions are necessary.

Advantages of the product

- Excellent material compatibility. Particularly suitable for use on sensitive surfaces made of plastic, acrylic glass or artificial leather.
- Comprehensive effective spectrum against bacteria, fungi and all viruses (virucidal with limited effectiveness against enveloped viruses, e.g. Vaccinia viruses, BVDV inclusive of HBV, HCV and HIV).
- VAH/DGHM registered. Tested in accordance with EN 14348, EN 14476.
- Dries quickly.

Ecology:

IQ B45 is easily biodegradable according to OECD guideline 301 D. The packing is made of polyethylene and can therefore be exploited both materially and thermally. To recycle, rinse the container with water.

Toxicology:

The acute oral and dermal toxicity (LD50) in rats is > 2000 mg/kg of body weight. Therefore, IQ B45 is classified as "non-toxic". In accordance with the results of an in vitro test, IQ B45 is classified as skin-irritating. If the preparation comes into contact with the mucous membranes of the eye, it must also be treated as an irritant. There is no risk of inhalation if the preparation is used as specified.

Microbiological Effectiveness

IQ B45 has the following characteristics of performance:

Application	Concentration	Time
Surface disinfection (according to VAH/DGHM)	2%	15 min.
Bacteria inclusive of MRSA ₁) and fungi ₂)	2%	15 min.
TB bacteria	2%	1 hour
Vaccinia viruses, BVDV incl. HBV, HCV and HIV ₃)	1%	1 hour

1) Effectiveness against bacteria implies effectiveness against MRSA.

2) Test for Candida albicans.

3) In accordance with RKI recommendation (Bundesgesundheitsbl. 47, 62 – 66, 2004).



**B45 Ready-to-use
Spray 500 mL**



**B45 Ready-to-use
Spray 1L**



B45 Ready-to-use 1L



**B45 Ready-to-use
Jerry Can 5L**



**B45 Ready-to-use
Stackable Jerry Can 10L**

Product Name	U/M ml	Pack Size	Art. No.
B45 Ready-to-use Spray Bottle	500 ml	06 pcs/ box	I020080016
B45 Ready-to-use Spray Bottle	1000 ml	09 pcs/ box	I020080119
B45 Concentrate Round Bottle with cap	1000 ml	06 pcs/ box	I020080006
B45 Concentrate Jerry can	5000 ml	04 pcs/ box	I020080054
B45 Concentrate Stackable Jerry Can	10000 ml	Single can	I020080010

Meets the Highest Standards of Hygiene

A safe disinfectant even if the concentrated tablet has been mentioned to be corrosive, once dissolved into the recommended solution TOSIL® is gentle to common surfaces and not irritating to skin. You can work safely with the solution. The Safety Data Sheet and product information sheets provide you with more informations and reports. In the recommended solution TOSIL® IS not harmful, cytotoxic or carcinogenic. It does not accumulate in the environment and is readily biodegradable.

TOSIL MED

TOSIL® med is a highly effective, safe and sustainable surface disinfectant - inline with EU Standards - and active against microorganisms. It is effective against a broad spectrum of viruses, bacteria, spores, fungi and yeast. TOSIL® med acts through the mechanism of oxidation and destroys microorganisms irreversibly. Thus it ensures that microbial resistance will not occur. Disinfecting with TOSIL® med helps to protect against infections.



Chemical description:
Sodium-p-toluenesulfonylchloramide (Chloramin-T)

TOSIL HITEC

TOSIL® hitec is a safe and sustainable surface disinfectant for all types of surfaces, equipments and machineries which can be rinsed, sprayed or wiped down. It's broad spectrum germicidal activity acts through the mechanism of oxidation and destroys microorganisms irreversibly. Thus it ensures that microbial resistance will not occur. Disinfection with TOSIL hitec is recommended in all fields where impeccable hygienic conditions are necessary.



Chemical description:
Sodium p-toluenesulfonylchloramide (2.5g each tablet)

	TOSIL® med	TOSIL® Hitec
Bacteria	MRSA, Enterococcus hirae, Pseudomonas aeruginosa, Escheria coli, Legionella pneumophila, Salmonella sp., Proteus vulgaris.	MRSA, Enterococcus hirae, Pseudomonas aeruginosa, Escheria coli, Legionella pneumophila, Salmonella sp., Proteus vulgaris.
Viruses	HIV, Norovirus, Alphavirus, H7N1, H1N1, Human Rotavirus, Adeno Virus, Hepatitis B Virus, Polio Virus, Corona Virus.	HIV, Norovirus, Alphavirus, H7N1, H1N1, Human Rotavirus, Adeno Virus, Hepatitis B Virus, Polio Virus.
Fungi	Candida albicans	Candida albicans
Spores	Clostridium difficile	Clostridium difficile

Product Name	U/M tab	Pack Size	Art. No.
Tosil Med	2.5 g/tab	25 pcs/ box	I010100025
Tosil Hitec	2.5 g/tab	25 pcs/ box	I010020025

TOSIL PLUS WIPES

Material friendly surface disinfectant combined with cleaning power!

RANGE OF APPLICATION

Acc. to MDD (Medical Device Directive): Disinfection and cleaning of alcohol-sensitive medical devices.

Acc. to BPD (Biocide Directive): Disinfection and cleaning of alcohol-sensitive medical inventory and all kinds of surfaces also in food industry and canteen kitchens.

APPLICATION

Wipe surfaces and devices until entirely moistened. After exposure time has elapsed wipe off surface with a single use cloth as required. TOSIL wipes plus are ready to use and can be extracted one by one from the dispenser by pulling them obliquely through the dispenser star.

For disinfection and cleaning

Alcohol-free wipes for disinfection and cleaning in hospital, primary healthcare and general practice. Used on all types of surfaces, medical devices and non-invasive appliances / instruments.

ENVIRONMENTAL INFORMATION

The products of Inoquest Labs L.L.C. are manufactured according to modern, safe and environmentally friendly processes. Due to the compliance with high quality standards the constant and very good product quality can be guaranteed.

PRODUCT PROPERTIES

Alcohol-free, ready-to-use disinfection wipes for disinfection and cleaning of medical inventory and all kinds of surfaces in medical areas as well as in all areas of the food industry and in canteens. Especially suitable for alcohol sensitive materials and equipments. TOSIL wipes plus active substance Sodium-ptolenesu fonchloramide is free from alcohol, aldehyde and phenols.

TOSIL wipes plus provide excellent cleaning results with no toxicity.

TOSIL
plus wipes



Tosil wipes Plus
100 sheets



Tosil wipes Plus
200 sheets

Product Name	U/M ml	Pack Size	Art. No.
Tosil Wipes Plus 100 Sheets	100 Sheets	20 pcs/ box	1070020012
Tosil Wipes Plus 200 Sheets	200 Sheets	12 pcs/ box	1070020020



INOQUEST LABS
Disinfectant and Hygiene Products



HAND DISPENSERS

Complementary accessories for disinfectants and soaps facilitate the day-to-day handling and dosage of the products used. Our dispenser range allow targeted application and dispensing of precise quantities, which optimise effectiveness and consumption.

INTRODUCTION

The COVID-19 pandemic revealed the importance of hand hygiene for infection control. Hand sanitizer dispensers are now commonplace in the entrance areas of office buildings and public places.

The demand for hand hygiene products has been unprecedented. We have invested heavily over the past year to meet those demands. To protect as many people as possible against the coronavirus and other infectious diseases, these investments will continue.

The World Health Organization (WHO) and the Robert Koch Institute in Germany recommend thorough and regular hand hygiene to strengthen infection protection as approximately 80 percent of all infections are transmitted by the hands. This fact underlines the importance of a proper hand hygiene infrastructure. On an average weekday we spend more than 8 hours outside our homes. We go to the work, grocery shopping, or to the hairdresser, we meet with friends in the evening, go to the cinema or have dinner with our partner. Frequently touched objects like door handles, rails, or chair backs can be a potential transmission source of harmful pathogens.

When this time of “freedom” comes again, we must be given opportunities to wash or sanitize our hands, particularly in high traffic areas and places where many people come together where germs can be transmitted quickly. Easily accessible sanitizer dispensers can offer effective protection against viruses and bacteria.



PRAESIDIO

The new standard for public
hand hygiene



PRAESIDIO

The PRAESIDIO® offers users a convenient opportunity for touchless hand hygiene to prevent infections. In addition to its elegant, high-quality design, the PRAESIDIO® offers a wide range of features: adjustable dosing quantity, hygienic spray chamber, integrated drip tray and viewing window for product filling level control. An optional stand enables flexible positioning of the PRAESIDIO® in corridors and highly frequented entrance areas for effective infection prevention.

The timeless, intuitive design of the PRAESIDIO® makes it the ideal candidate for hand hygiene in high-traffic public areas. This robust dispenser was designed to support effective infection prevention thanks to its touchless operation and user-friendly experience. The patented spray nozzle, spray chamber, integrated drip tray, and fill-level window complete the PRAESIDIO®. Optional stands ensure easy, safe access in highly frequented areas.



KX SERIES

**RELIABLE
AUTOMATIC
SENSING**



KX Touch Free

- ADA Compliant
 - Touchless operation
 - Low-profile
- Easy one-handed bottle exchange
- Locking or non-locking cover
- Compatible with all 38 mm X Series cartridge pumps
- LED indicator light, signals operation and low battery
- Efficient energy consumption: 60,000 – 90,000 activations
- High-grade ABS plastic



KX Manual

Sleek and compact, the KX series while it fits into small spaces. Compatible with all OPHARDT Hygiene 38 mm pumps, the KX Series provides industry-leading touchless and manual dispensing for liquid, foam, and spray.



- ADA Compliant
 - Low-force presser
 - Low-profile
- Easy one-handed bottle exchange
- Fits anywhere, saves wall space
- Easy one-handed bottle exchange
- Locking or non-locking cover
- Compatible with all 38 mm X Series pumps
- High-grade ABS plastic

ERGONOMIC AND EASY TO USE



NEPTUNE SERIES

Elegance Meets Flexibility

This elegantly styled dispenser series offers a wide range of dispensing options and an unmatched



Neptune Touch Free Dispenser

- ADA Compliant
 - Touchless operation
 - Low-profile
- Easy one handed bottle exchange
- Locking or non-locking cover
- Compatible with all 38 mm X Series Cartridge pumps
- LED indicator light signals operation and low battery
- Efficient energy consumption: 60,000 – 90,000 activations



Sensor



Manual

Neptune Manual Dispenser

- ADA Compliant
 - Low-force presser
 - Low-profile
- One-handed dispensing
- Easy one handed bottle exchange
- Locking or non-locking cover
- Compatible with all 38 mm X Series Cartridge pumps

D10I MANUAL DISPENSER



With the introduction of the ingo-man® plus dispenser system in 2006, the flagship series from IQ was elevated. Being able to quickly remove the pump from the front of the dispenser significantly reduced maintenance time and effort.

In addition, a flexible silicone membrane now seals the different bottlenecks so that the product is always protected from contamination. Rapid pumping and optimized dosing are other advantages the plus series has over the ingo-man® classic series.



Manual

D10I Dispenser

- White plastic housing
- Plastic control lever
- Aluminium back panel, lockable
- For euro bottles with external diameter of the neck thread of 22 - 60 mm
- Plastic pump for soaps, disinfectants and lotions.
- For 1000 ml standard Euro bottle, can also be filled freely

POUCH DISPENSERS

Economic Triangle Pouch Dispensing System.

- White plastic with grey push button
- Sensor operated
- 4 X 1,5V battery Size C
- Easy to refill
- Lockable
- 1,000 ml



Sensor



Manual

Economic Pouch Dispensing System.

- White plastic with grey push button
- Manual operated
- Easy to refill
- Lockable
- 1,000 ml



IQ Stainless Steel

Surgical Scrub Sink & Transfer Trolley

► IQ Surgical Scrub Sink - Single Station



Art. No. I050200001

- Made of SS304G stainless steel with satin finish and most optimal RA value.
- Automatic hands free / elbow operated water and automatic / foot operated hand scrub.
- One water tap and one hand scrub tap.
- Curved joints as per HTM standards ensuring perfect hygiene and ease of cleaning.
- Openable back splash for easy refilling of hand scrub and routine maintenance.
- The unit is wall mounted for easy cleaning of the floor.

Dimension : 33"L x 24.5"W x 31.25"H

Water gun with 8 different tips (Optional)

ISO 2015-9001 Certified manufacture.



Scrub Tap



Scrub Inside

► IQ Surgical Scrub Sink - Two Station



Art. No. I050200002

- Made of SS304G stainless steel with satin finish and most optimal RA value.
- Automatic hands free / elbow operated water and automatic / foot operated hand scrub.
- Two water tap and one hand scrub tap.
- Curved joints as per HTM standards ensuring perfect hygiene and ease of cleaning.
- Openable back splash for easy refilling of hand scrub and routine maintenance.
- The unit is wall mounted for easy cleaning of the floor.

Dimension : 63"L x 24.5"W x 31.25"H



Scrub Tap



Scrub Inside

► IQ Surgical Scrub Sink - Three Station



Art. No. I050200003

- Made of SS304G stainless steel with satin finish and most optimal RA value.
- Automatic hands free / elbow operated water and automatic / foot operated hand scrub.
- Three water tap and two hand scrub tap.
- Curved joints as per HTM standards ensuring perfect hygiene and ease of cleaning.
- Openable back splash for easy refilling of hand scrub and routine maintenance.
- The unit is wall mounted for easy cleaning of the floor.

Dimension : 93"L x 24.5"W x 31.25"H



Scrub Tap



Scrub Inside

► IQ Basket Trolley with Eight Baskets



Wire net Basket with/
without silicon inserts



SPRI Baskets

Art. No. I050020001



- Manufactured from 304G stainless steel
- SS stopper for preventing basket fall during transportation
- Suitable for transportation of sterilization baskets
- TENTE wheels with diagonal brakes

■ Available in :

SINGLE(4 baskets) /

DOUBLE (8 baskets) column trolley

■ Dimension :

48"L x 17.5"B x 79"H (Double Column)

24.75"L x 17.5"B x 79"H (Single Column)

- Suitable for sterilization devices
- Easy stack / stock design

Basket Type 22.75 - 1" L x 15.5" W x 8"H

Basket Type 22.75 - 2" L x 15.5" W x



► IQ Closed Transportation Trolley



Art. No. I050010003

- Manufactured from 304G stainless steel
- Lockable door and double side pushing handle
- Adjustable horizontal shelves - 3 nos
- TENTE wheels with diagonal brakes
- Visual indicator for Sterile / processed and Contaminated/ not processed material (Optional)
- Flexible production in different designs and dimensions

Dimension : 48"L x 24"W x 50"H



IQ Waste Bins



Ingo-man[®] Waste Bins

ingo-man®

WASTE BINS

Made of stainless steel. available with Anti-Fingerprint Coating, optimized for use in medical facilities, rollable, with foot-operated opening mechanism, with internal waste bag holder, self-closing, soft-close lid, recessed lid to prevent hand contact (minimizes contamination), fire-inhibiting, free-standing with wheels.



NWB 24

▶	24 Liter	Art. No.
	NWB 24 E	Stainless steel 1416587
	NWB 24 E AFP	Stainless steel with AFP Coating 1418560
▶	50 Liter	Art. No.
	NWB 50 E	Stainless steel 1416588
	NWB 50 E AFP	Stainless steel with AFP Coating 1419160
▶	70 Liter	Art. No.
	NWB 70 E	Stainless steel 1416589
	NWB 70 E AFP	Stainless steel with AFP Coating 1419161


*AFP - Anti-Fingerprint Coating



NWB 50



NWB 70

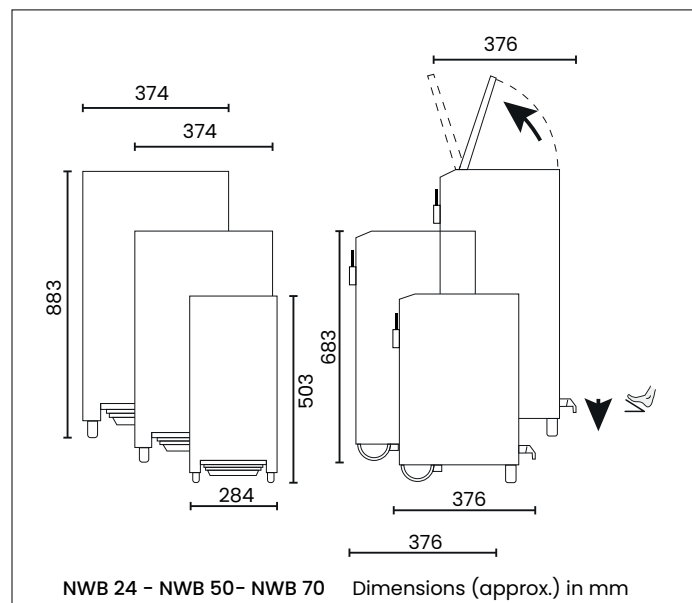



Rectangular Waste Bin

Stainless Steel Waste Bin




RWB 30L

RWB 20L

RWB 12L

- Damper system to control lid closing speed
- Automatic soft close
- Quiet closing
- With foot pedal
- With inner bucket
- Stainless steel
- Warranty up to one year

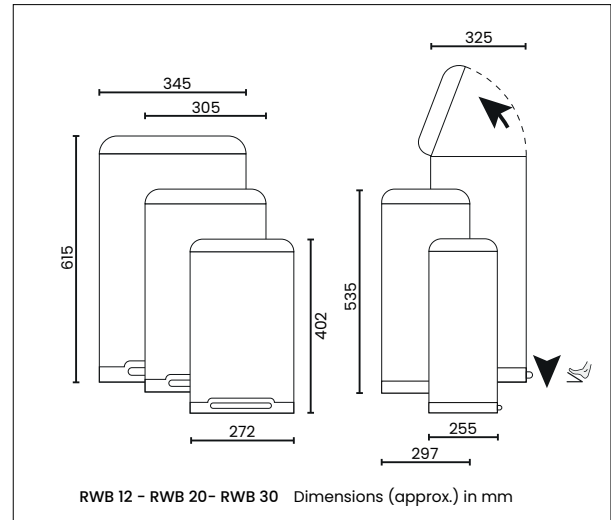

Powder Coated Steel Waste Bin

- Damper system to control lid closing speed
- Automatic soft close
- Quiet closing
- With foot pedal
- With inner bucket
- Powder coated steel
- Warranty up to one year



Anti-Fingerprint



Slowly & Quietly



The unique air holes,
extraction of garbage bags easy

Powder coated, with lid and foot pedal Waste bin made of steel with powder coating. For general and medical waste.

Product	Material	Dimension	Art. No.
IQWBS 12L	Stainless steel	272×255×402mm	IQWBS06SS12L
IQWBS 20L	Stainless steel	305×297×535mm	IQWBS06SS20L
IQWBS 30L	Stainless steel	345×325×615mm	IQWBS06SS30L
IQWBS 12L Black	Powder coated steel	272×255×402mm	IQWBS06SB12L
IQWBS 30L Black	Powder coated steel	345×325×615mm	IQWBS06SB30L
IQWBS 12L Yellow	Powder coated steel	272×255×402mm	IQWBS06SY12L
IQWBS 30L Yellow	Powder coated steel	345×325×615mm	IQWBS06SY30L
IQWBS 12L Red	Powder coated steel	272×255×402mm	IQWBS06SR12L
IQWBS 30L Red	Powder coated steel	345×325×615mm	IQWBS06SR30L
IQWBS 30L White	Powder coated steel	345×325×615mm	IQWBS06SW30L



Rectangular Waste Bin 50 L

Stainless Steel Waste Bin 50L



RWB 50L



With foot pedal




With inner bucket

Waste bin made of brushed stainless steel. For general and medical waste.

Product	Material	Dimension	Art. No.
IQWBS 50L	Stainless Steel	475×390×645mm	IQWBS07SS50L

- Damper system to control lid closing speed
- Automatic soft close
- Quiet closing
- With foot pedal
- With inner bucket
- With transport handles
- Stainless steel

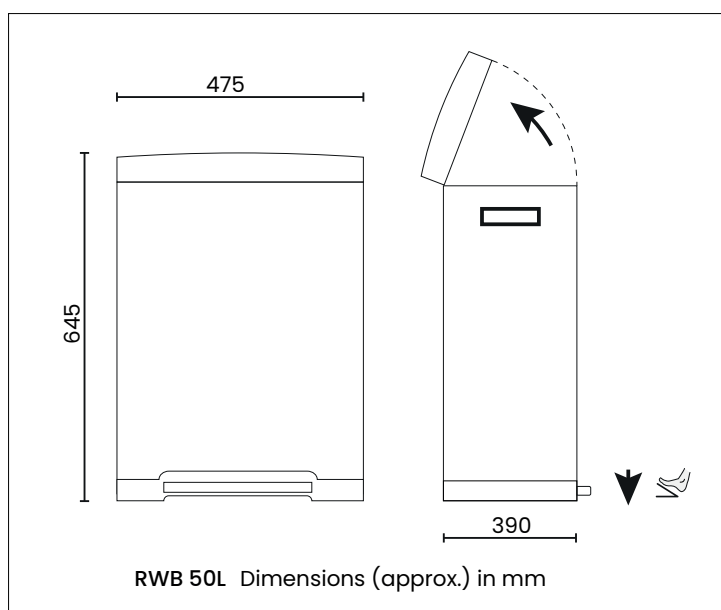



Round Waste Bin

Round Staineless Steel Waste Bin




RWB 12L

RWB 20L

Product	Material	Dimension	Art. No.
IQWBR 12L	Stainless steel	315x 250x 390mm	IQWBR07SS12L
IQWBR 20L	Stainless steel	360x 293x 430mm	IQWBR07SS20L

- With foot pedal
- With inner bucket
- With transport handles
- Stainless steel


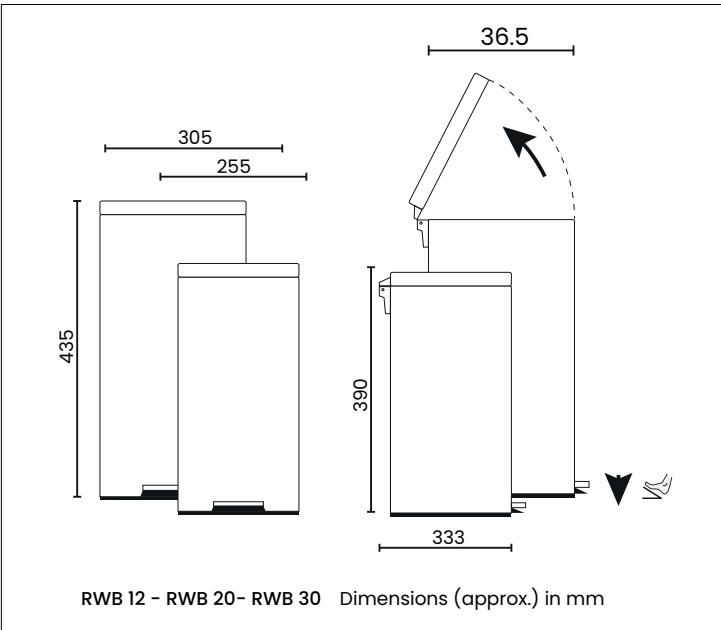




IQ Paper Towel Dispenser

01

HTU 1 Paper Towel Dispenser Touchless

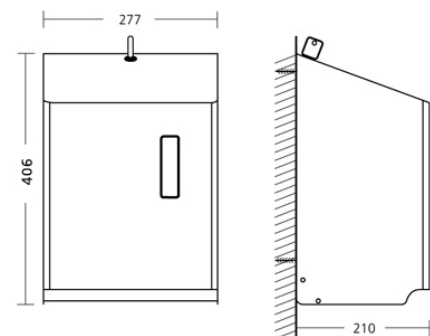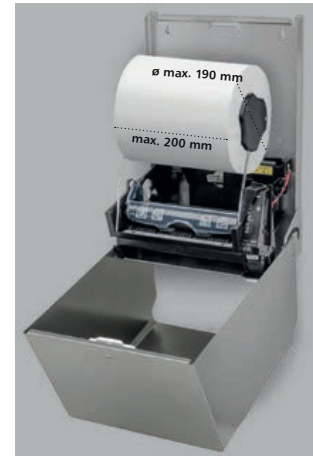

HTU | dimensions (approx.) in mm

Housing made of stainless steel, lockable. with window for monitoring filling level

Paper feed: Touchless paper feed Paper remains completely within the dispenser until a hand is detected underneath the removal point - delay and run-on length adjustable

Filling Material: paper rolls
max. roll width 200 mm, max. roll diameter 190 mm,
min. inside diameter for paper withdrawal 38 mm

Power Supply: with 4 D-batteries for a total of 6 Volts

Dimensions: W 278x H 406 x D 211 mm

HTU 1 E AFP


Stainless steel with Anti-Fingerprint Coating

Art. No.


3400351

02


Z fold / Interfold Paper Dispenser

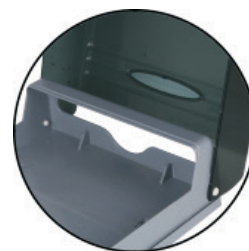

Easy to take

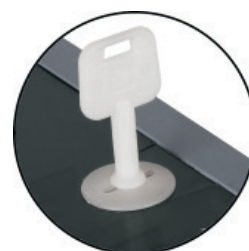

Key

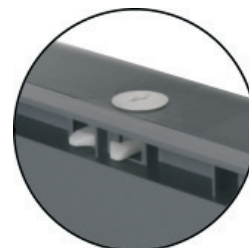


Lockable

Art. No.: I040130010

1. Easy to install
2. Lockable design, cover to prevent vandalism/theft/tampering of contents

Material	Product Size (MM)	PCS/CTN	Carton Size (CM)	N.W./G.W.(KG)
ABS	375*95*270	10	58.5*40*58.5	14/12.5

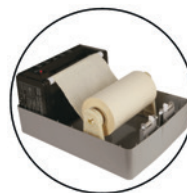

Easy to fit paper

Inside movement
for AutocutAdjustable Length and 2 type
working way: sensor or automatic


Size of paper

Art. No. : I040340050


Sensor towel paper dispenser Advantage:

- 1, Support DC and battery (4x1.5V / D size)
- 2, Long battery life (more than 10000times)
3. Sheet length can be adjustable : S/M/L (20/25/30)
4. 2 working type: sensor (put hand for sensor then the paper comes out) or automatic (pull the paper then come out, always have paper no need sensor)
5. Sensor type: if the paper is not took off, will stop working to save paper, need pull the paper, then sensor working
- 6, Suitable for paper: core size more than 3.7cm, length: 20cm, diameter: 20cm
7. Normally roll paper can use around 700 times (25cm)
8. Security key and lock
9. Color Options : White , Translucent Black

Material	Product Size (MM)	PCS/CTN	Carton Size (CM)	N.W./G.W.(KG)
ABS	233*315*318	2	26*40*66	10/8

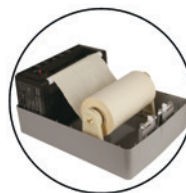

Easy to fit paper

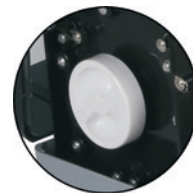

Button is available

Inside movement
for Autocut

Size of paper


Art. No. : I040340010

Autocut towel paper dispenser Advantage:

- 1.No need power, just pull the paper by hand
- 2.Long time use (more than 300000times)
3. Suitable for paper: core size more than 3.7cm, length:20Cm, diameter: 20cm
4. Normally 1 roll paper can use around 700 times (25cm)
5. Security key and lock
6. Color Options : White , Translucent Black

Material	Product Size (MM)	PCS/CTN	Carton Size (CM)	N.W./G.W.(KG)
ABS	233*315*318	2	26*40*66	8.5/7.5

Our Quality Certifications

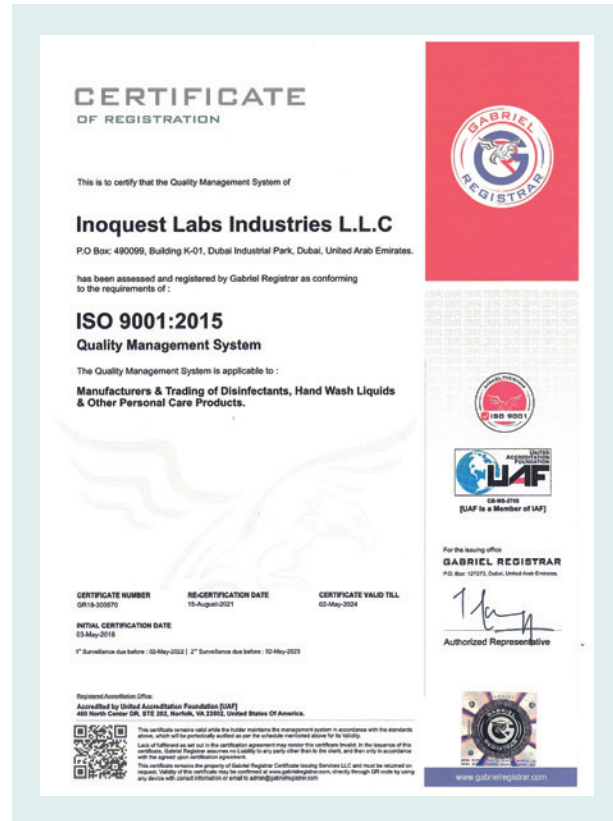

We are a recognized product manufacturer who follows GMP and ISO criteria.

This means that our products are subjected to extremely strict quality controls and regulations.

We have a long history of upholding this quality for our products and clientele. From R&D through quality control, we've got you covered. Quality inspections are extremely important to us, especially during the manufacturing process.

We maintain consistent standards by continuous refining of our technique and comprehensive product testing. Our products are also evaluated by external institutes and specialists for their efficacy and skin tolerance.

Manufactured & Distributed By:



DIN ISO 9001:2015 Certified Company
GMP & DIN ISO 22716:2007

***Your product — our expertise
Are our products just right for you?***

Drop us a line or give us a call!

+971 4 450 8696

Products & Services

- Hospital & Hygiene Products
- Homecare Products
- Disinfectants & Detergents
- Disinfection Equipments
- Industrial Hygiene Products
- Private Label / OEM Products

Address

Inoquest Labs Industries L.L.C.
Dubai Industrial City, WH K-01,
P.O. Box 49 00 99, Dubai, UAE.

Phone

(+971) 4 450 8696

Online

Mail : info@inoquestlabs.com

Website : www.inoquestlabs.com